KINYOUN COLOR KIT

(Tan-Thiam-Hok and Devulder method)

ITALIANO

DESCRIZIONE

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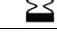
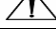
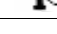
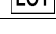
ENGLISH

DESCRIPTION

KINYOUN COLOR KIT (Tan-Thiam-Hok and Devulder method) is a kit for colouring alcohol-acid resistant bacteria in clinical samples or culture media. The Kinyoun colouring (Tan-Thiam-Hok and Devulder method) is an alternative to colouring alcohol-acid resistant bacteria that is simpler to administer than hot Ziehl Neelsen colouring. Colouring can be carried out using the following methods:

- Tan-Thiam-Hok,
- Devulder, hot version
- Devulder, cold version

CONTENT OF THE PACKAGES

The reagents are contained in plastic bottles that are heat sealed and provided with drip caps. Each kit contains:

- 1 bottle containing 250 ml of Kinyoun Solution
- 1 bottle containing 250 ml of Gabett Solution
- 1 instruction sheet

PRINCIPLE OF THE METHOD

Kinyoun colouring (Tan-Thiam-Hok and Devulder method) is based on the ability of phenol to enable basic fuchsia to penetrate inside the bacterial cell. Certain bacteria, one they have been coloured with basic fuchsia, if they are treated with alcohol-acid decolouring solution remain coloured red (alcohol-acid resistant). Others lose their purple colour and take on the colour of the methylene blue contrast (bacteria sensitive to decolouring with alcohol-acid).

COMPOSITION

KINYOUN COLOR KIT		
Content		
Kinyoun solution	Basic fuchsia	40 g
	Ethyl alcohol 95°C	200 ml
	Phenol	80 ml
	Distilled water	1000 ml
Gabett solution	Methylene blue	10 g
	Sulphuric acid	200 ml
	Ethyl alcohol	300 ml
	Distilled water	500 ml

COLLECTION OF SAMPLES

The samples to be subjected to colouring of the acid-resistant alcohols mainly consist of clinical material and microbial cultures.

TEST PROCEDURE

Preparation and fixing

Use clean plates. Swab the culture or pathological material. Dry in the air and fix by passing rapidly over the flame. Other fixing methods can be used.

Staining

During all phases of the above methods, the preparations must be washed in water free of alcohol-acid resistant species.

TAN-THIAM-HOK METHOD

- Cover the plate with Kinyoun solution and leave in contact for 3 minutes.
- Wash with water for 30 seconds.
- Cover the plate with Gabett solution and leave in contact for 1 minute.
- Wash with water and dry.
- Observe the preparation under the microscope using the immersion lens.

HOTDEVULDERMETHOD

- Cover the plate with Kinyoun solution and heat for 10 minutes.
- Wash with water.
- Cover the plate with Gabett solution. Leave in contact for 5 minutes.
- Wash with water and dry.
- Observe the preparation under the microscope using the immersion lens.

COLD DEVULDER METHOD

- Immerse the plate in the Kinyoun solution for about 3 hours.
- Wash with water.
- Cover the plate with Gabett solution. Leave in contact for 5 minutes.
- Wash with water and dry.
- Observe the preparation under the microscope using the immersion lens.

INTERPRETATION OF THE RESULTS

Mycobacteria and other bacteria that are resistant to alcohol and acid show up red whilst background material and other bacteria that are sensitive to colouring with acid-alcohol show up blue.

QUALITY CONTROL

Each batch of KINYOUN COLOR KIT (Tan-Thiam-Hok and Devulder method) is subjected to quality control according to the three methods described using an ATCC 14468 *Mycobacterium smegmatis* culture to test for bacteria that are resistant to alcohol and acid (show up red) and an ATCC 25923 *Staphylococcus aureus* culture to test for bacteria that are discoloured by alcohol-acid (show up blue).

LIMITS AND WARNINGS

The intensity of the colouring of the bacteria that are resistant to alcohol and acid is less intense than with traditional Ziehl Neelsen colouring but the granulation is more clearly visible.

PRECAUTIONS

The KINYOUN COLOR KIT (Tan-Thiam-Hok and Devulder method) contains substances that are classified as hazardous by current legislation: consult the safety datasheet before use. KINYOUN COLOR KIT (Tan-Thiam-Hok and Devulder method) is a kit for bacterial colouring that is only for diagnostic *in vitro* professional use. It must be used in the laboratory by properly trained operators using approved aseptic and safety methods for pathogenic agents.

STORAGE

Store KINYOUN COLOR KIT at 10-25°C in its original packaging. Do not store it near sources of heat and do not expose it to excessive temperature variations. If properly stored, the KINYOUN COLOR KIT (Tan-Thiam-Hok and Devulder method) can be used until the expiry date shown on the label. Do not use after this date. Dispose of it if there are signs of deterioration (solutions change colour or there are signs of large precipitates).

DISPOSAL OF USED MATERIAL

After use, the coloured plates with the KINYOUN COLOR KIT (Tan-Thiam-Hok and Devulder method) and the material that has come into contact with the sample must be decontaminated and disposed of in accordance with current laboratory techniques for the decontamination and disposal of potentially infected material.

BIBLIOGRAPHY

- B.G. Metchock, F.S. Nolte & R.J. Wallace, 1999. *Mycobacterium*, p.399-437 In P.R. Murray, E.J. Baron, M.A. Pfaller, F.C. Tenover and R.H. Tenover (eds.), *Manual of Clinical Microbiology*, 7th ed. American Society for Microbiology, Washington, D.C.
- Devulder B. Ann. De l'Institut Pasteur de Lille. 1963, **14**, 229-231.
- Gabett H.S. Lancet 1987. **1**, 757.
- Kinyoun J.J. Amer. J. Public Health, 1915, **5**, 867.
- Murray, P.R. (ed.) 1999. *Manual of Clinical Microbiology*, 7th ed. American Society of Microbiology, Washington, D.C.
- Tan-Thiam-Hok Am. Rev. Resp. Dis., 1962, **85**, 753.

PRESENTATION

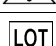
Product	REF	
KINYOUN COLOR KIT (Tan-Thiam-Hok and Devulder method)	80380	2 x 250 ml

TABLE OF SYMBOLS

SYMBOL	MEANINGS	SYMBOL	MEANINGS
	In Vitro Diagnostic Medical Device		Do not reuse
	Manufacturer		Contains sufficient for <n> tests
	Catalogue number		Fragile, handle with care
	Use by		Caution, consult accompanying documents
	Temperature limitation		Batch code

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